beurer

FC 76

DE Mikrodermabrasionsgerät Gebrauchsanweisung	2	TR Mikrodermabrazyon cihazı Kullanım kılavuzu	51
EN Microdermabrasion device Instructions for use	12	RU Прибор по уходу за кожей лица Инструкция по применению.....	61
FR Appareil de microdermabrasion Mode d'emploi.....	21	PL Urządzenie do mikrodermabrazji Instrukcja obsługi.....	72
ES Aparato de microdermoabrasión Instrucciones de uso.....	31		
IT Apparecchio per microdermoabrasione Istruzioni per l'uso.....	41		

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89231 Neu-Ulm

Germany

Read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

Table of contents

1. Getting to know your instrument	13
2. Signs and symbols	14
3. Intended use	14
4. Warnings and safety notes	15
5. Unit description	16
6. Use	16
6.1 Carrying out a skin test	17
6.2 Using the device.....	17
6.3 Application areas.....	18
7. Maintenance and cleaning.....	19
8. Replacement parts and wearing parts	19
9. Disposal.....	20
10. Technical data.....	20
11. Warranty/service	20

Dear customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for applications in the areas of heat, weight, blood pressure, body temperature, pulse, gentle therapy, massage, beauty, baby and air.

Please read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

With kind regards,
Your Beurer team

Included in delivery

- 1x FC 76
- 1x precision attachment
- 1x coarse attachment
- 1x fine attachment
- 20x single-use filter
- 1x mains adapter
- These instructions for use

WARNING

- The device is only intended for domestic/private use, not for commercial use.
- This device is not suitable for persons under the age of 18.
- This device may be used by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely and are fully aware of the consequent risks of use.
- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children unless supervised.
- If the mains connection cable of this device is damaged, it must be disposed of. If it cannot be removed, the device must be disposed of.

1. Getting to know your instrument

The FC 76 enables professional facial peels with microdermabrasion technology. A high-quality sapphire coating is used to gently remove excess skin cells and the state-of-the-art vacuum technology noticeably improves skin texture and enhances circulation. The treatment also simulates cell renewal for a finer and smoother complexion.

2. Signs and symbols

The following symbols appear in these instructions for use and on the device:

	WARNING Warning instruction indicating a risk of injury or damage to health
	IMPORTANT Safety note indicating possible damage to the device/accessory
	DANGER This device must not be used near to water or in water (e.g. in a wash basin, or in the shower or bath) – risk of electric shock!
	Note Note on important information
	Observe the instructions for use.
	The device is double-protected.
	Disposal in accordance with the Waste Electrical and Electronic Equipment EC Directive – WEEE
	Manufacturer
	Only use indoors.
	Dispose of packaging in an environmentally friendly manner

3. Intended use

The device is only to be used on the face (apart from the eye area and mouth) of adults. Do not use the device on animals!

The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage resulting from improper or careless use.

4. Warnings and safety notes

WARNING

- Danger of suffocation by packaging material. Keep packaging material away from children.
- Before use, inspect the device and accessories for visible damage. If you have any doubts, do not use the device and contact your retailer or the specified Customer Services address.
- Keep the mains adapter away from water.

DANGER

This device must not be used near to water or in water (e.g. in a wash basin, or in the shower or bath) – risk of electric shock!

WARNING

Do **NOT** use the device

- on cracked skin
- on open wounds
- in the eye area or on the mouth
- on any part of the body other than the face (except for eye area and mouth)
- if you have a skin disease (e.g. herpes, vascular diseases or skin growths)
- if you have skin cancer
- if you are suffering from skin irritation
- immediately after sunbathing or if you have sunburn
- if you are taking medication that contains steroids
- if you are taking blood thinners
- if you have multiple sclerosis, diabetes or an autoimmune disorder

WARNING

For this reason, the device must only be operated as follows only with the supplied power supply unit at the mains voltage indicated on the power supply unit.


IMPORTANT

Only specialist personnel may perform repairs on electrical devices. Improper repairs may subject users to considerable danger. In the event of a fault or damage, have the device repaired by a certified specialist workshop.

WARNING

- Only use the device in dry rooms (e.g. never in the bath, sauna).
- Never reach for a device that has been submerged in water. Disconnect the mains adapter immediately.
- Do not drop the device or subject it to any impacts.
- Do not clamp or pinch the mains cable.
- Ensure that the mains cable and mains adapter do not come into contact with water, steam or other liquids.

5. Unit description

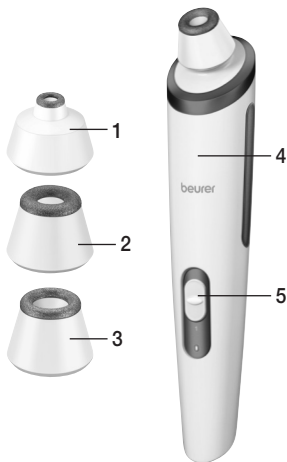1. Precision attachment with fine sapphire coating, suitable for precise treatment of small areas (e.g. in nose area)
2. Fine attachment with fine sapphire coating, suitable for sensitive skin
3. Coarse attachment with coarse sapphire coating, suitable for intensive treatment of problem areas
4. Hand-held device
5. ON/OFF switch and adjustment of intensity levels 1-2

6. Use

Note

Before you start using the device, make sure that your face is clean, dry and free of any traces of products such as powder, cream or make-up.

6.1 Carrying out a skin test

Carry out a skin test 24 hours before using the device. Test the device at the lowest intensity level on small, inconspicuous areas of your neck. If after 24 hours the tested area is red or painful, do not use the device.

WARNING

- Do not use the device more than once a week.
- Immediately after use, the treated skin may be slightly reddened and in very rare cases there may be small spots of bleeding under the skin.
- For one day before and one day after use, avoid activities that cause you to perspire, such as saunas, tanning beds and sun exposure.


Note


The device can also be used twice a week for non-sensitive problem skin (with a break of at least 48 hours between use).

6.2 Using the device


1. Connect the mains adapter to the device and plug the adapter into the socket.

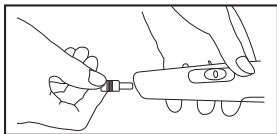
2. Insert a single-use filter into the head of the device.

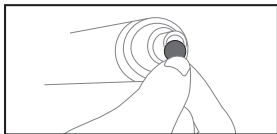3. Place the attachment you wish to use on the tip of the hand-held device.

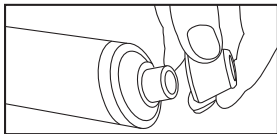4. To switch on the device, move the ON/OFF switch to the desired intensity level.

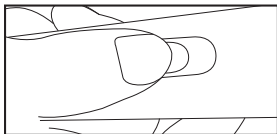
5. Place the hand-held device on your face and move it over the skin. Make sure that you do not move the hand-held device over one part of the face more than once. Further details on treatment can be found in Section "6.3 Application areas".

i Note

The pressure you apply when using the device should be comfortable and not painful. Some slight reddening should be expected after use.

6. When you have completed the application, switch off the device using the ON/OFF switch and disconnect the mains adapter from the socket..

7. Change the single-use filter after each use.
To do this, pull the attachment off from the top. Use tweezers to remove the single-use filter.

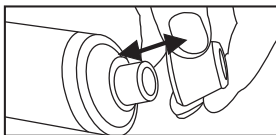
6.3 Application areas

i Note

Ensure that the part of your face you are treating is taut when you move the hand-held device over the skin. You can do this by pulling the skin taut to the side or downwards with your free hand.

Forehead:

Move the hand-held device slowly and evenly outwards from the centre of the forehead. Make sure that you do not treat any part of your skin more than once.

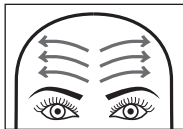Temples:

Move the hand-held device outwards towards your hairline. Repeat this to the left and right of the nose. Ensure that the device is at least 1 cm away from your eyes.

Nose:

Move the hand-held device downwards from the centre at the top of the bridge of your nose. Repeat this to the left and right of the nose.

Cheeks:

Move the hand-held device outwards across your cheeks.

Around the mouth:

Move the hand-held device along your upper lip from the centre to the left and right. Move the hand-held device upwards from your bottom lip to the left and right. Make sure that you do not touch your lips.

7. Maintenance and cleaning

- Never submerge the device in water or other liquids.
- Disinfect the attachments with a standard disinfectant after each use.
- Change the single-use filter after each use.
- Use a dry cloth to clean the device of any dirt. Do not use any abrasive cleaning products.
- Protect the device from knocks, damp, dust, chemicals, marked temperature fluctuations, and nearby sources of heat (ovens, heaters).

8. Replacement parts and wearing parts

Replacement parts and wearing parts are available from the corresponding service address (according to the service address list).

Please state the corresponding order number.

Replacement part	Item number and/or order number
10 single-use filters	163.660
Precision attachment	163.767
Coarse attachment	163.766
Large attachment	163.765

9. Disposal

For environmental reasons, do not dispose of the device in the household waste at the end of its useful life. Dispose of the unit at a suitable local collection or recycling point.

Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.

10. Technical data

Power unit voltage	Input: 100–240V AC / 50/60 Hz Output: 5V DC / 1A
Dimensions	190 x 40 x 39 mm
Weight	~128 g

11. Warranty/service

In case of a claim under the warranty please contact your local dealer or the local representation which is mentioned in the list “service international”.

In case of returning the unit please add a copy of your receipt and a short report of the defect.

The following warranty terms shall apply:

1. The warranty period for BEURER products is either 3 years or- if longer- the country specific warranty period from date of Purchase.
In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Repairs (complete unit or parts of the unit) do not extend the warranty period.
3. The warranty shall not be valid for damages because of
 - a. improper treatment, e.g. nonobservance of the user instructions.
 - b. repairs or tampering by the customer or unauthorised third parties.
 - c. transport from the manufacturer to the consumer or during transport to the service centre.
 - d. The warranty shall not be valid for accessories which are subject to normal wear and tear (cuff, batteries etc.).
4. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

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