beurer

HK 55 Easyfix

DE Multifunktionsheizkissen Gebrauchsanweisung	2	IT Termoforo multifunzione Istruzioni per l'uso.....	29
EN Multifunctional heat pad Instruction for Use.....	9	TR Çok işlevli ısı yastığı Kullanma Talimatı	36
FR Coussin chauffant multifonctions Mode d'emploi	15	RU Электрогрелка Инструкция по применению	42
ES Almohadilla eléctrica multifunción Instrucciones para el uso	22	PL Wielofunkcyjna poduszka elektryczna Instrukcja obsługi	50

CE

Inhalt

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Germany

Contents

<p>1. Items included in the package10</p> <p style="padding-left: 20px;">1.1 Description10</p> <p>2. Important instructions</p> <p style="padding-left: 20px;">Retain for future use10</p> <p>3. Intended use12</p> <p>4. Operation12</p> <p style="padding-left: 20px;">4.1 Safety12</p> <p style="padding-left: 20px;">4.2 Initial use12</p> <p style="padding-left: 20px;">4.3 Switching on12</p>	<p>4.4 Setting the temperature12</p> <p>4.5 Automatic switch-off13</p> <p>4.6 Switching off13</p> <p>5. Cleaning and maintenance13</p> <p>6. Storage14</p> <p>7. Disposal14</p> <p>8. What if there are problems?14</p> <p>9. Technical data14</p> <p>10. Warranty/service14</p>
--	--

Explanation of symbols

	Read the instructions!		Machine-washable on an extra-delicate cycle at 30°C
			Do not bleach
	Do not insert pins!		Do not tumble dry
			Do not iron
	Do not use folded or rucked!		Do not dry-clean
	Not to be used by very young children (0 - 3 years).	<div style="border: 1px solid black; padding: 5px; width: fit-content; margin: 0 auto;"> <p style="font-size: 8px; margin: 0;">OEKO-TEX® CONFIDENCE IN TEXTILES STANDARD 100 Special Article 06.0.43510 HOHENSTEIN HTTI <small>Tested for harmful substances, www.oeko-tex.com/standard100</small></p> </div>	
 	Dispose of packaging in an environmentally friendly manner		
	WARNING: Warning of risks of injury or health hazards		


	<p>CAUTION: Safety information about possible damage to appliance/accessories.</p>
	<p>NOTE: Important information.</p>

1. Items included in the package

- 1 Heating pad
- 1 Controller
- 1 Extension strap
- 1 Instructions for use

1.1 Description

- 1. Mains plug
- 2. Mains cable
- 3. Controller
- 4. Illuminated display for temperature settings
- 5. Sliding switch for ON/OFF and temperature settings
- 6. Plug-in coupling

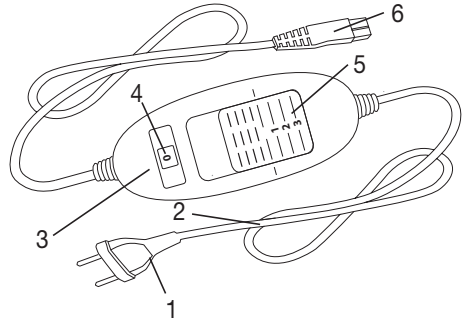2. Important instructions

Retain for future use

WARNING

- Non-observance of the following notes may result in personal injury or material damage (electric shock, skin burns, fire). The following safety and hazard information is not only intended to protect your health and the health of others, it should also protect the product. For this reason, pay attention to these safety notes and include these instructions when handing over the product to others.
- This heating pad must not be used by people who are not sensitive to heat or by vulnerable people who may not be able to react to overheating (e.g. diabetics, persons with skin changes related to illness or persons with scarred skin in the area where the pad is to be applied, after taking any pain-relieving medication or alcohol).
- This heating pad may not be used by very young children (0-3 years old) as they are unable to respond to overheating.
- This heating pad may not be used by young children (3-8 years old) unless the controller has been set by a parent or responsible adult and the child has been given adequate instructions on how to use the heating pad safely.

- This heating pad may be used by children over the age of 8 and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised and have been instructed on how to use the heating pad safely, and are fully aware of the consequent risks of use.
- Children must not play with the heating pad.
- Cleaning and user maintenance must not be performed by children unless supervised.
- This heating pad is not designed for use in hospitals.
- Do not insert needles.
- Do not use when folded or bunched up.
- Do not use if wet.
- This heating pad may only be used in conjunction with the controller specified on the label.
- The electrical and magnetic fields emitted by this heating pad may interfere with the function of a pacemaker. However, they are still well within the limits: electrical field strength: max. 5000 V/m, magnetic field strength: max. 80 A/m, magnetic flux density: max. 0.1 millitesla. Please consult your doctor and the manufacturer of your pacemaker before using this heating pad.
- Do not pull, twist or make sharp bends in the cables.
- Please check this heating pad frequently for signs of wear and tear or damage.
If any such signs are evident, if the heating pad has been used incorrectly or if it no longer heats up, it must be checked by the manufacturer before being switched on again.
- If the mains connection cable of this device is damaged, it must be disposed of. If it cannot be removed, the device must be disposed of.
- While the heating pad is switched on, the following conditions must be observed:
 - Objects must not be placed on the pad (e.g. suitcase or washing basket).
 - Heat sources, such as hot water bottles, heating pads or similar, must not be placed on the pad.
- The electronic components in the controller warm up when the heating pad is in use. For this reason, the controller must never be covered or placed on the heating pad when it is in use.
- Always observe the notes on operation (section 4), cleaning and maintenance (section 5) and storage (section 6).
- If you should have any questions about using our devices, please contact Customer Services.

3. Intended use

CAUTION

This heating pad is only designed to warm up the human body.

This heating pad is suitable for universal use. For example, on the stomach, back, joints or neck. The heating pad can be rolled for use on the neck area.

The extendible hook-and-loop fasteners make the device very easy to attach to the body. Thanks to the rubber strap the heating pad always fits perfectly.

4. Operation

4.1 Safety

CAUTION

The heating pad is fitted with a SAFETY SYSTEM. This sensor technology provides protection against overheating across the entire surface of the heating pad with an automatic switch-off in the event of a fault. If the SAFETY SYSTEM has switched off the heating pad, the temperature settings are no longer illuminated when switched on. Please note that for safety reasons, the heating pad can no longer be operated after a fault has occurred and must be sent to the specified service address.


4.2 Initial use

CAUTION

Ensure that the heating pad will not bunch up or become folded during use.

- First, connect the controller to the heating pad by inserting the plug-in coupling.
- Then plug the mains plug into the socket.

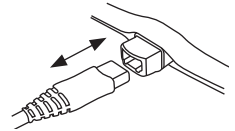4.3 Switching on

To switch the heating pad on, set the sliding switch for ON/OFF and temperature settings to setting 1, 2 or 3. The temperature settings (4) are illuminated when switched on.

4.4 Setting the temperature

Level 0: OFF

Level 1: minimum heat

Level 2: medium heat

Level 3: maximum heat

NOTE

The fastest way to warm up the heating pad is to initially set the highest temperature setting.

NOTE

These heating pads have a fast heating function, which allows the pad to warm up faster in the first 10 minutes.

WARNING

If the heating pad is being used over several hours, we recommend that you set the lowest temperature setting on the controller to avoid overheating the heated body part, which may lead to burns to the skin.

4.5 Automatic switch-off

This heating pad is equipped with an automatic switch-off function. This turns off the heat supply approx. 90 minutes after the initial use of the heating pad. The temperature setting display then starts to flash. So that the heating pad can be switched back on, the sliding switch for ON/OFF and temperature settings (5) must first be set to setting "0" (OFF). After about 5 seconds it is possible to switch it on again.

4.6 Switching off

To switch the heating pad off, set the sliding switch for ON/OFF and temperature settings to setting "0" (OFF). The temperature settings are then no longer illuminated.

NOTE:

If the heating pad is not in use, switch the sliding switch for ON/OFF and temperature settings to setting "0" (OFF) and unplug the mains plug from the socket. Then disconnect the controller from the heating pad by unplugging the plug-in coupling.

5. Cleaning and maintenance

WARNING

Before cleaning, always remove the mains plug from the socket first. Then disconnect the controller from the heating pad by unplugging the plug-in coupling. Otherwise there is a risk of electric shock.

CAUTION

- The controller must never come into contact with water or other liquids, as this may cause damage.

To clean the controller, use a dry, lint-free cloth. Do not use any chemical or abrasive cleaning agents.

CAUTION

- Please note that the heating pad may not be chemically cleaned, wrung out, tumble dried, put through a mangle or ironed. Otherwise the heating pad may be damaged.

This heating pad is machine-washable. Set the washing machine to a particularly gentle wash cycle at 30°C (wool cycle). Use a delicate laundry detergent and measure it out according to the manufacturer's instructions.

CAUTION

- Please note that frequent washing of the heating pad has a negative effect on the product. The heating pad should therefore be washed in a washing machine a maximum of 10 times during its life.

Immediately after washing, reshape the heating pad to its original dimensions while it is still damp and spread it out flat on a clothes horse to dry.

CAUTION

- Do not use pegs or similar items to attach the heating pad to the clothes horse. Otherwise the heating pad may be damaged.
- Do not reconnect the controller to the heating pad until the plug-in connection and the heating pad are completely dry. Otherwise the heating pad may be damaged.

WARNING

- Never switch the heating pad on to dry it! Otherwise there is a risk of electric shock.

6. Storage

If you do not plan to use the heating pad for a long period of time, we recommend that you store it in the original packaging. For this purpose, disconnect the controller from the heating pad by unplugging the plug-in coupling.

CAUTION

Please allow the heating pad to cool down before storing. Otherwise the heating pad may be damaged. To avoid sharp folds in the heating pad, do not place any objects on top of it while it is being stored.

7. Disposal

For environmental reasons, do not dispose of the device in the household waste at the end of its useful life. Dispose of the unit at a suitable local collection or recycling point. Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.

8. What if there are problems?

Problem	Cause	Solution
The temperature settings are not illuminated while - the controller is properly connected to the heating pad - the mains plug is connected to a working socket - the sliding switch is set to 1, 2 or 3	The safety system has switched the heating pad off permanently.	Send the heating pad and the controller for servicing.

9. Technical data

See the rating label on the heating pad.

10. Warranty/service

In case of a claim under the warranty please contact your local dealer or the local representation which is mentioned in the list “service international”.

In case of returning the unit please add a copy of your receipt and a short report of the defect.

The following warranty terms shall apply:

1. The warranty period for BEURER products is either 3 years or- if longer- the country specific warranty period from date of Purchase.
In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Repairs (complete unit or parts of the unit) do not extend the warranty period.
3. The warranty shall not be valid for damages because of
 - a. improper treatment, e.g. nonobservance of the user instructions.
 - b. repairs or tampering by the customer or unauthorised third parties.
 - c. transport from the manufacturer to the consumer or during transport to the service centre.
 - d. The warranty shall not be valid for accessories which are subject to normal wear and tear (cuff, batteries etc.).
4. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

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